



<b>DECLARATION OF KOON-WAH LEONG UNDER 37 C.F.R. §1.131</b>  Address to: Commissioner for Patents P.O. Box 1450 Alexandria, VA 22313-1450	Application Number	09/593,827
	Confirmation Number	7505
	Filing Date	June 13, 2000
	First Named Inventor	Sherry X. Guo
	Examiner	R. Gitomer
	Group Art	1623
	Attorney Docket No.	LIFE008/LFS-101

**RECEIVED**  
SEP 25 2003

TECH CENTER 1600/2900

This Declaration and the attached Exhibit are being submitted in conjunction with the Applicants' Response to the Office Action dated June 16, 2003.


I, Koon-wah Leong, Ph.D., do hereby declare as follows.

1. I am an inventor of the invention claimed in the above captioned application, along with Sherry Guo.
2. I have been asked to declare and provide factual evidence in support of a reduction to practice of a storage stable composition containing nylon and 10-(carboxymethylaminocarbonyl)-3,7-bis(dimethylamino)phenothiazine before October 26, 1999.
3. As evidenced by Exhibit A, we reduced to practice a storage stable composition containing nylon and 10-(carboxymethylaminocarbonyl)-3,7-bis(dimethylamino)phenothiazine prior to October 26, 1999. The dates have been redacted from Exhibit A. All redacted dates are prior to October 26, 1999.

4. Exhibit A consists of a photocopy of three pages from my lab notebook, showing an experiment that I performed called "Trial 3". The following is written at the end of page 651-105: "Rgt II was coated on OneTouch Nylon membrane at 65°C for 10 min". As is apparent from the middle of page 651-105, Rgt II contains B dip. B dip contains component DA-67, which corresponds to 10-(carboxymethylaminocarbonyl)-3,7-bis(dimethylamino)phenothiazine. As such, this page discloses a nylon membrane coated with 10-(carboxymethylaminocarbonyl)-3,7-bis(dimethylamino)phenothiazine. Pages 651-106 and 651-107 show test results, and concludes with "Conclusion: very good calibration curve was obtained". According to my lab notebook, the membrane was coated with DA-67 prior to October 26, 1999. The date of the signature at the bottom of the page is also prior to October 26, 1999.
5. The evidence provided in Exhibit A establishes that we produced a storage stable composition containing nylon and 10-(carboxymethylaminocarbonyl)-3,7-bis(dimethylamino)phenothiazine prior to October 26, 1999.
6. I hereby declare that all statements made herein of my own knowledge are true and that all statements made on information and belief are believed to be true; and further that these statements were made with the knowledge that willful false statements and the like so made are punishable by fine or imprisonment, or both, under Section 1001 of Title 18 of the United States Code and that such willful false statements may jeopardize the validity of the application or any patents issued thereon.

Respectfully submitted,

Date: September 9, 2003

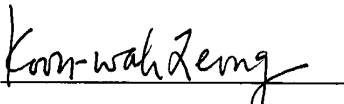  
Koon-wah Leong, Ph.D.

Attachments: Exhibit A

LIFESCAN INCORPORATED

EXHIBIT A 1/3

No


651-105 P E

DATE PROJECT

09/593,827

SUBJECT

From Page No. 104

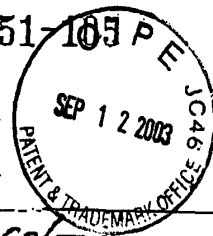

When the above strips were tested, no blue color was observed due to the inhibition of KAO by Bis-MAPS-C3.

Trial 3 :

RECEIVED

SEP 25 2003

Rgt I : 200 mg/ml of protease TR

4% mannitol

0.5% Triton X-100

5mm Ca(Rc)<sub>2</sub>

in 0.1M EPPS, pH = 8.0.

TECH CENTER 1600/2

Rgt I was coated on Rubicon strip version 3. 2-layer coating was applied. (coated on

Rgt II : A dip of 200 u/ml KAO

1mg/ml HRP

4% mannitol

1% PVP 360K

in 0.1M phosphate, pH = 7.5.

B dip of 1mm DA-67 in 70% methanol, 30% H<sub>2</sub>O

Rgt II was coated on One-touch Nylon membrane @ 65°C for 10 min.

Yellow soln. was observed when B dip was coated.

One-touch Nylon membrane was mounted on Sure-strip membrane strip holder.

Data On Disk No. Cont'd. on page

Signed Date

Read and Understood Symon Qian Date

No. 651-106

4089423825

T-682 P.004/005 F-185


LUMESCAN INCORPORATED

DATE \_\_\_\_\_ PROJECT \_\_\_\_\_

EXHIBIT A 2/3 09/593,827

SUBJECT \_\_\_\_\_

From Page No. 104

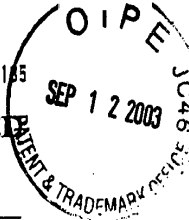

Testing Procedures:



- ① 100ul of control serum was added to Rubicon strips coated w/ Rgt I. The strips were then put on 37°C heating plate for 20 min
- ② 25ul digested serum pooled for 5 Rubicon strips were then added to ~~one~~ strips coated w/ Rgt I.
- ③ The color product was monitored on ~~plate~~

RECEIVED

Results

SEP 25 2003

TECH CENTER 1600/2900

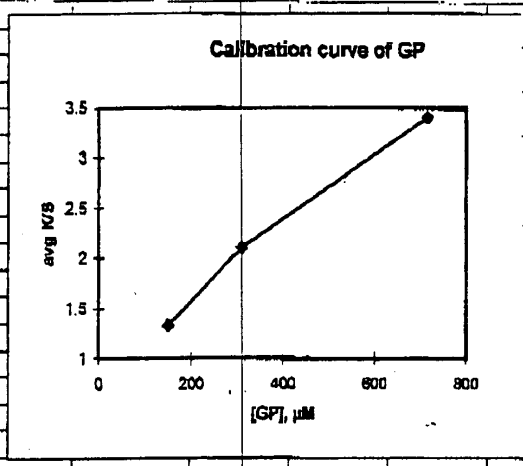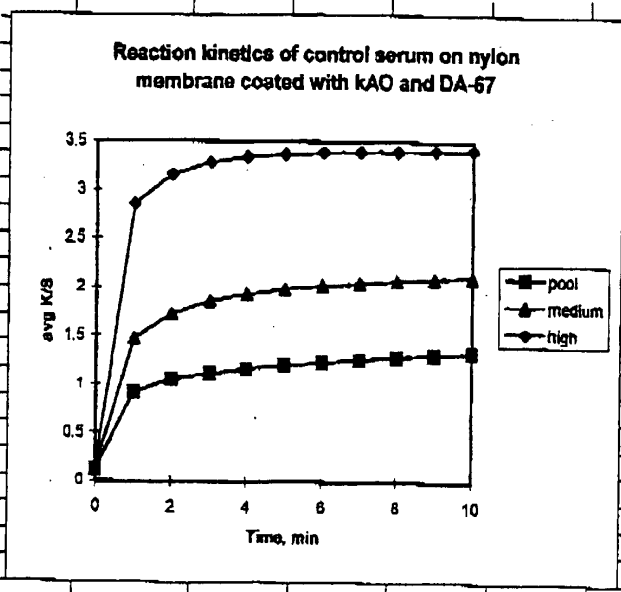

Data On Disk No. \_\_\_\_\_ Cont'd. on page \_\_\_\_\_

Signed \_\_\_\_\_ Date \_\_\_\_\_

Read and Understood Suyue Qian Date \_\_\_\_\_

LIFESCAN INCORPORATED

EXHIBIT A 3/3

Nº


651-107

DATE \_\_\_\_\_ PROJECT \_\_\_\_\_

09/5/93, 827

SUBJECT \_\_\_\_\_

From Page No. 104

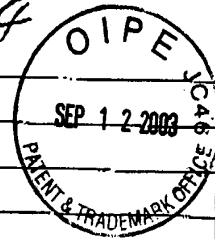

		λ=660nm					
	[Time] min	R%1	R%2	K/S1	K/S2	avg K/S	
Standard		87.95					
Pool serum 070797	Time, min					pool	
150uM by GivPro	0	62	62	0.116452	0.116452	0.116452	
	1		28.12		0.918694	0.918694	
	2	26.54	25.47	1.016649	1.090444	1.053548	
	3	25.55	24.81	1.084697	1.139366	1.112032	
	4	24.85	24.18	1.136322	1.188725	1.162524	
	5	24.33	23.71	1.176726	1.227365	1.202045	
	6	23.94	23.33	1.208255	1.259813	1.234034	
	7	23.59	23.03	1.237492	1.286231	1.261862	
	8	23.3	22.74	1.262423	1.312469	1.287448	
	9	23.05	22.5	1.284447	1.334722	1.309585	
	10	22.81	22.25	1.306071	1.358441	1.332258	
Medium 1 090397	Time, min					medium	
310 uM by GivPro	0	62	62	0.116452	0.116452	0.116452	
	1	20.2	22.22	1.578248	1.381325	1.468788	
	2	18.38	19.63	1.812248	1.645272	1.72878	
	3	17.67	18.45	1.918005	1.802277	1.860141	
	4	17.31	17.82	1.975054	1.894938	1.934995	
	5	17.05	17.45	2.017801	1.96258	1.98519	
	6	16.87	17.19	2.048191	1.994618	2.021404	
	7	16.72	17.02	2.074031	2.02282	2.048425	
	8	16.59	16.88	2.096814	2.046485	2.07185	
	9	16.48	16.78	2.116381	2.063638	2.090009	
	10	16.38	16.89	2.134403	2.079256	2.106829	
High 090397	Time, min					high	
716 uM by GivPro	0	62	62	0.116452	0.116452	0.116452	
	1	12.47	14.03	3.071973	2.633942	2.852957	
	2	11.78	12.84	3.303382	3.018896	3.161139	
	3	11.53	12.16	3.394163	3.172642	3.283403	
	4	11.37	12	3.454387	3.226667	3.340527	
	5	11.33	11.88	3.469713	3.268154	3.368933	
	6	11.28	11.83	3.489024	3.285893	3.387358	
	7	11.26	11.81	3.496797	3.29275	3.394774	
	8	11.25	11.79	3.500694	3.299832	3.400263	
	9	11.23	11.77	3.50851	3.308938	3.407724	
	10	11.22	11.76	3.512428	3.310501	3.411484	

RECEIVED

SEP 25 2003

TECH CENTER 1600/200

Conclusion:

Very good calibration curve was obtained.  
However, the interference still need to be checked.

Data On Disk No. \_\_\_\_\_ Cont'd. on page \_\_\_\_\_

Signed \_\_\_\_\_ Date \_\_\_\_\_

Read and Understood Ernie Qian Date \_\_\_\_\_